


## PERFCONT2E Perfusion Module

The PERFCONT2E perfusion module is an option for the SYSTEM2/10EH *Wave Bioreactor*. The BASE2/10EH is placed on the PERFCONT2E platform. This provides a continuous weight measurement. Tare functions allow the monitoring and control of the culture volume. The BASE2/10EH keypad and LCD display allow various modes to be programmed to feed or remove preset amounts of media or harvest, and for automatic perfusion control. The current flowrate, the amount added, and the cumulative amounts added/removed are displayed.

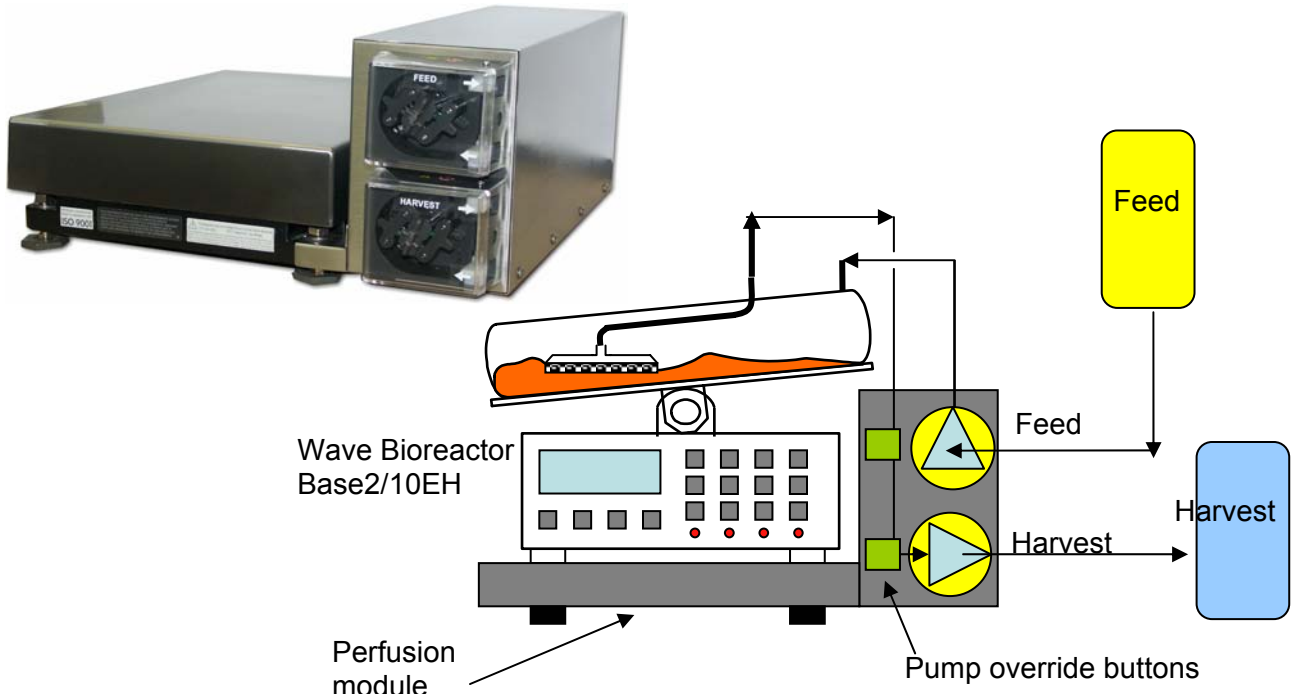

<b>Performance:</b>	Stainless-steel easy-to-clean housing. Two internal peristaltic pumps for tubing from 0.5mm to 4.8mm diameter. Load cell weight sensor 0 to 15Kg (Accuracy: 0.02% F.S. Max 55Kg.) Control and power interface for BASE2/10EH. Pump override buttons for test and priming.
<b>Dimensions:</b>	350mm x 330mm x 150mm (13.75 x 12.65 x 6.05 inch)
<b>Weight:</b>	10Kg (22lb)
<b>Power:</b>	Interface cable to BASE2/10EH.
<b>Options:</b>	None